# BEAUTIFLY

SlimSense

Body

series


Waga analityczna  
Figure Scanner

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[f /beautifly.eu](https://www.facebook.com/beautifly.eu) [@beautifly.eu](https://www.instagram.com/beautifly.eu) [g beautifly.eu](https://www.beautifly.eu)

[/beautiflyeu](https://www.youtube.com/channel/UC...) [@beautiflyeu](https://www.tiktok.com/@beautiflyeu)


f1

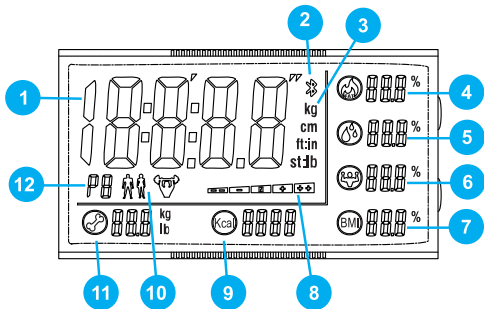f2

f3

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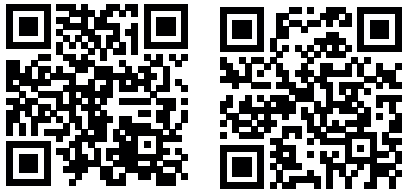- [www.beautify.eu](http://www.beautify.eu)

- [www.youtube.com/beautifyeu](http://www.youtube.com/beautifyeu)

Type of scale:	Electronic
Surface type:	ITO glass
Maximum weight:	180kg
Power supply:	4xAAA batteries
Features:	measurement memory, body analysis: weight, body fat content, body water, muscle measurement, bone mass, BMI, Kcal

#### Safety

Before using the product, read the instructions and follow the recommendations contained therein during the use of the device.

#### Display Description [page 2 - f1]

1. Display data / body weight
2. Bluetooth
3. Units
4. Body fat percentage (%)
5. Body water (%)
6. Muscles (%)
7. BMI
8. Body fat level: low/standard/moderate/overweight/obesity
9. Kcal
10. Male/Female/Athlete
11. Bone mass
12. User profile

#### How to use

1. Before first use, install batteries. To do this, open the cover on the back of the scale and insert the batteries according to the polarity indication (+/-).
2. Place the scale on a hard and flat surface, DO NOT use it on a carpet or soft surface.
3. Make sure that the feet and the surface of the scale are dry to avoid slipping.
4. Stand straight with bare feet, covering the surface of the scale evenly [page 2 - f2]
5. Press the "SET", "▲" or "▼" button on the scale to turn it on, and the scale will go to the user parameter settings. The device stores personal parameters for 9 people (P1-P9).
6. Select your profile.

Note: the scale has a memory function, there is no need to set user data each time. If the weight difference between you and other users is within +/-3 kg during one weighing, the scale will identify all users and display the corresponding user number.

7. Stand still for a few seconds and do not get off the scale until the data appears on the LCD display. The device shows 7 indicators: weight, body fat content, body water, muscles, bone mass, BMI, Kcal.

Note: The scale will automatically turn off after 20 seconds if no other operation is performed.

When the scale is not connected to Bluetooth - its symbol on the LCD screen will flash, and the LCD display will show the user number from the range "P1" - "P9"; When the scale is connected to Bluetooth, its symbol will be displayed continuously, and the user number "P0" will appear on the LCD display.

#### Setting personal parameters:

1. Press <SET> to turn on the scale.
2. Press ▲ or ▼ to toggle within P1-P9 and confirm by pressing <SET>
3. To set personal parameters press ▲ or ▼ and select one of the options, then confirm by pressing <SET>: male - female - athlete / height / age
4. Press ▲ or ▼ to select the unit of mass (kg/lb/st/lb). Press <SET> to confirm.

#### Hazards and warnings:

1. Do not use replacement parts that are not approved by the manufacturer.
2. Non-professional service should not disassemble and replace parts of the device without permission.
3. If the scale operates incorrectly, check the batteries. If necessary, replace the batteries.
4. If the percentage of body fat is less than 5% or more than 50% the scale shows an error - check if the user information is correct.
5. The scale is a precision measuring device. Never jump or hit the scale.
6. Store away from water, heat, and extreme cold.
7. Make sure the device is not exposed to direct sunlight or placed near heaters and other heat sources.
8. The data provided on the scale is for reference only. Before starting any diet or exercise program, consult a doctor.
9. The scale is not a toy
10. The device is not suitable for pregnant women and people with a pacemaker

#### Cleaning and maintenance:

Before cleaning, turn off the product. To clean the dirty surface, use a soft cloth soaked in alcohol or glass cleaner. Do not clean the product and accessories with corrosive detergents, such as benzene, thinners. Store in a cool, dry place when not in use. Also avoid exposure to direct sunlight or a humid environment, otherwise it may cause color fading or damage to the product.

#### More information about the company and products [page 2 - f3]

- [www.beautify.eu](http://www.beautify.eu)
- [www.youtube.com/beautifyeu](https://www.youtube.com/beautifyeu)

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after 13.08.2005.

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[EN] The manufacturer declares under its sole responsibility that the SlimSense product is in compliance with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, Directive 2014/53 of the European Parliament and of the Council /EU of 16 April 2014 on the harmonization of the laws of the Member States relating to the making available on the market of radio equipment and repealing Directive 1999/5/EC, and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances. The full content of the declaration of conformity is available at: [www.beautifyfy.eu/deklaracje](http://www.beautifyfy.eu/deklaracje)

**Importer:** Beautify Sp. z o.o. Plac Bankowy 2, 00-095 Warszawa, Polska  
**Model:** SlimSense