DIFORO
FAIRY 2
ProSeries

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Vobe sp. z o. o. declares that the device Diforo is in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following internet address www.diforo.pl in the declarations tab.

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EN-The crossed-out wheeled bin symbol indicates that in the European Union after the end of the product's useful life it must be disposed of in a separate, dedicated collection point. Do not dispose of these products together with unsorted general household waste

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DIFORO
FAIRY 2
ProSeries

TERMS OF WARRANTY

1. This warranty is provided by VOBÉ sp. z o.o., ul. Rudna Mała 167A, 36-054 Mrowla, registered by the District Court for the Capital City of Warsaw in Warsaw, 12th Commercial Division of the National Court Register under the KRS number: 0000800189, with the assigned numbers NIP: 525-279-92-70, REGON: 384166183, hereinafter referred to as the Guarantor.
2. The warranty period is 24 months and is counted from the date of delivery of the product to the original Buyer. The warranty does not cover the accessories and the battery.
3. This warranty, to the extent and within the period specified therein, provides the Buyer with free of charge removal of any defects in the advertising product. The liability under the guarantee covers only defects resulting from the reasons inherent in the sold item.
4. Product defects will be removed by the Guarantor's Authorized Service, hereinafter referred to as the Service.
5. In order to use the rights resulting from this warranty (filing a complaint), the Buyer is obliged to fill in the complaint form or contact in another way with the Guarantor or the Service Center.
6. After making the notification, the Buyer is obliged to deliver the advertised product to the Website. We recommend that you include complete (factory) equipment with accessories with the product, as it may facilitate and speed up the warranty procedure.
7. The buyer is obliged to attach to the advertised product this warranty card together with the proof of purchase, as well as a description of the suspected product defects.
8. We recommend keeping the original packaging of the product in order to protect the product against damage that may occur during transport to the Service Center.
9. The guarantor will respond to the submitted complaint within 14 calendar days from the date of receipt by the Service Center of the product being complained about, along with the required set of documents referred to in point 7 of this guarantee.
10. The basic method of settling the complaint is to repair the product, restoring its utility value. If: a) the defect cannot be removed, b) the performance of the warranty repair is not possible within 14 calendar days from the date of receipt of the product by the Service Center, c) during the warranty period, three repairs qualified as warranty repairs have been made, but the product still has defects. The buyer is entitled to replace the product with a new or equivalent product, free from defects, of at least the same class. The above does not exclude a different method of settling the complaint, in a manner agreed jointly by the Guarantor and the Buyer.
11. In the event of non-compliance of the complaint with the facts or the inability to repair the product, the Guarantor and the Buyer shall agree on the manner and date of settling the complaint.
12. If the complaint is found to be unjustified, the Guarantor will return the product to the Buyer together with a written statement of the reasons for refusing to accept the complaint.
13. To the extent permitted by the provisions of applicable law, this warranty does not cover: a) the interaction of the Product with hardware or applications not supplied by the Guarantor, b) defects caused by the Buyer's use of the product contrary to its intended use or information contained in the product manual, c) damage the product resulting from random events, not resulting from the reasons inherent in the sold item, such as fire, flood, network overload, incorrect network voltage, lightning, warfare, etc. d) use of the product in conditions of significant pollution, higher humidity, extreme thermal conditions or environmental, or sudden changes to such conditions, e) defects and damage to the product resulting from the use of devices or applications not supplied by the Guarantor, such as chargers, fasteners, etc. f) opening, repairing, as well as changes made to the product by persons other than those

representatives of the Website.

14. To the extent permitted by applicable law, this warranty will not be enforceable if: a) the warranty card has been altered on its own,

b) any warranty seals have been removed from the product on its own accord.

15. To the extent permitted by applicable law, the Guarantor is not responsible for the loss of software or data stored in the device during repair, we recommend that you create a backup copy before returning the product for repair, as these data may be lost during the repair.

16. If, in the performance of his obligations, the Guarantor has provided the Buyer with an item free of defects instead of a defective item or has made significant repairs to the item covered by the warranty, the warranty period starts anew from the moment of delivery of the item free of defects or return of the repaired item. If the Guarantor replaced part of the item, the above provision shall apply to the replaced part. In other cases, the warranty period is extended by the time during which the Buyer could not use it due to a defect in the item covered by the warranty.

17. This warranty does not exclude, limit or suspend the rights of the Purchaser resulting from the non-compliance of the goods with the contract, as well as other rights resulting from the provisions of applicable law.

EN- WARRANTY / CLAIMS The product comes with a 24-month manufacturer's warranty. The warranty does not cover defects or damage caused by the fault of the user, or as a result of improper use (or misuse). The guarantee does not cover mechanical damage and caused as a result of the use of accessories other than those included in the kit. The guarantee does not cover consumable items such as batteries, rechargeable batteries. The warranty on batteries is 6 months. Complaints should be submitted in writing, e-mail address: get@diforo.pl

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EN- Declaration of Conformity with regard to the Directive RE 2014/53/UE, RE 2014/53/UE, RoHS 2011/65/UE

Adres serwisu / service address

Vobe sp. z o.o., Rudna Mała, 167A, 36-054 Mrowla
mail: get@diforo.pl, mob. +48 453 271 100
www.diforo.pl

Including

Charging cable

Instruction manual

Features

Microcurrent

The microcurrent **allows** the **electrical stimulation** to penetrate deep into the subcutaneous tissue to the **muscles**, helping to repair the skin's **elastic fibers** and the **gelatinous layer** tissue, thereby stretching and reducing wrinkles.

Warm Head Surface

The role of warming is to **effectively relieve** eye fatigue and promote **blood circulation**. Adherence to use can reduce dark circles, bags under the eyes, edema, etc., and can promote the absorption of eye cream.

Red Light Function

Red light therapy, known as photobiomodulation (PBM), encompasses one of the latest methods of treatment and regeneration of the body, using the natural properties of light with a wavelength of 620-630nm. It contributes to accelerating regenerative and healing processes in the body. The result is faster repair of damaged tissues. The application of red light therapy is wide-ranging and includes the treatment of problems such as skin problems including acne and eczema, as well as promoting wound and scar healing. Red light is used in anti-ageing therapy as it stimulates collagen and elastin production, reducing wrinkles and improving skin elasticity.

Blue Light Function

LED blue light therapy, with a wavelength of 400-495nm, focuses on seborrhoea, which is essential for the skin and hair to function properly. Blue light is able to focus on seborrhoea and make it less active. It eliminates acne-causing bacteria, which supports the treatment of more serious types of acne. It helps to remove scars and supports anti-inflammatory effects.


Vibration function

This product uses an **all metal** body and oxidation process technology, and the **overall** texture is very good. With an **electronic smart sensor** chip, when the product body is held, the head touches the skin, the machine starts to work, and stops when it leaves. Human computer interaction design, easy to operate.

Smart touch switch

This product uses an **all metal** body and oxidation process technology, and the **overall** texture is very good. With an **electronic smart sensor** chip, when the product body is held, the head touches the skin, the machine starts to work, and stops when it leaves. Human computer interaction design, easy to operate.

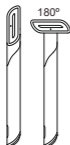
180 degree rotating head

Innovative 90 degree rotation **allows** the machine to better adapt to different parts of the eye.

EN

How to use?

It is recommended to use after cleaning the face for better results.

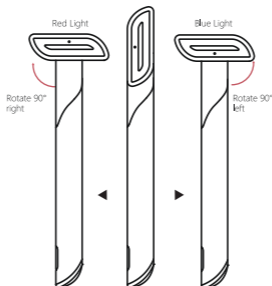- ① The head of the product rotates in two directions, each rotation can be made by 90°. During contact with the skin, Fairy will start to produce EMS currents and vibrations. When the treatment is complete, the head must be turned in the opposite direction to that before starting to switch off the device.
- ② This product uses electronic smart touch chip, intelligent switch. When holding the product body, the head touches the skin machine and starts working, leaving to stop working.
- ③ Except the eye area, this product is also applicable to the face and neck. Using as below method.

Eye area

Face

Neck

About charging

Battery Reminder

- When the battery is low, the indicator flashes quickly
- When charging, the indicator light is in the state of breathing light
- When fully charged, the indicator light is always on

Charging time is 120 minutes, continuous use time is 90 minutes.

Charging method

Plug one end of the charging cable into the charging port on the bottom of the product and the other end to connect the charger or the USB port of the computer. The indicator light flashes and charges.

 No more than 12 hours for charging time.

NOTES

- Do not disassemble or modify this product.
- Please follow the instructions to rotate the head.
- Do not rinse directly with water.
- After use, wipe with a rag or a wet tissue.
- Avoid direct sunlight when not in use for a long time.
- Keep this product out of the reach of children.

Product
Parameters:

Producer	DIFORO
Model	Fairy2
Material	Aluminium+PC
Size	177x16x19 mm
Net weight	28g
Battery capacity	300mAh
Charging time	120min
Working time	30min
Input voltage	5V
Current	0.3A
Power	1,5W
Certificate	CE RoHS

Importer:

Vobe sp. z o.o.
Rudna Mała 167A
36-054 Mrowla
NIP: 5252799270

www.diforo.pl

Made in China

EN